


THE A.P.C. SERIES

MODEL ASC275 & ASC375

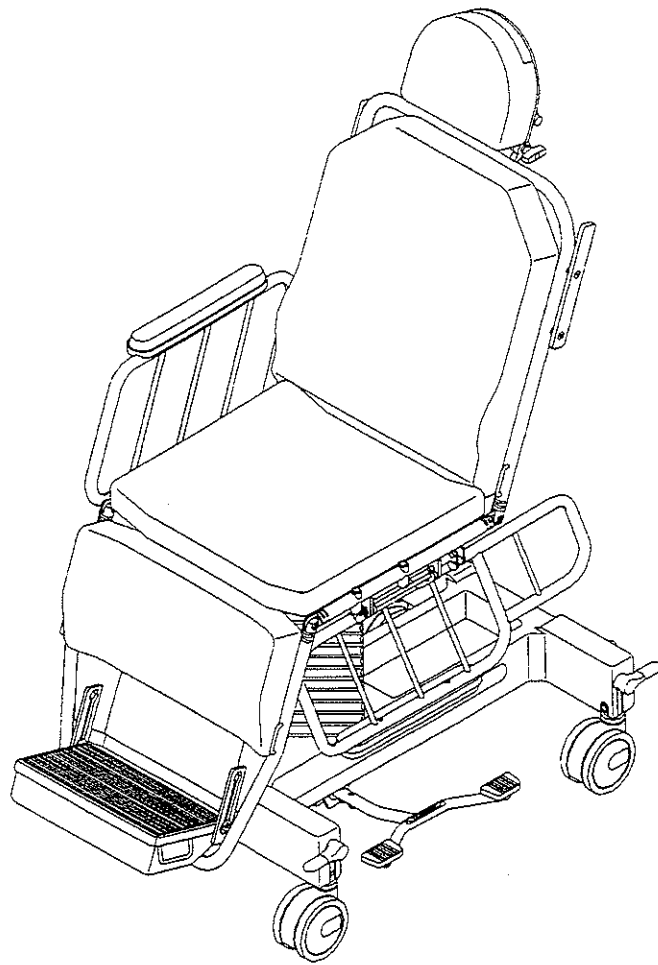

OPERATING AND SERVICE MANUAL

927 LAKE ROAD P.O. BOX 710 MEDINA, OHIO 44258
(216) 723-3271 1-800-428-7833

P/N 069071-65 (9604)



TABLE OF CONTENTS

PAGE 1	- UNPACKING INSTRUCTIONS
PAGE 2	- GENERAL INSTRUCTIONS - CHAIR SPECIFICATIONS
PAGE 3	- BRAKE AND STEER OPERATION - TOP ADJUSTMENT
PAGE 4	- TOP ADJUSTMENT - CONTINUED - USING THE SURGICAL LEGS
PAGE 5	- INDEPENDENT LEG SECTION OPERATION - ADJUSTABLE FOOTREST OPERATION
PAGE 6	- HINGING RAIL OPERATION - OPERATING THE A.S.C. HEADREST
PAGE 7	- BASE ASSEMBLY
PAGE 8	- SEAT SECTION ASSEMBLY
PAGE 9	- FOOT SECTION ASSEMBLY
PAGE 10	- BACK SECTION ASSEMBLY - MODEL 275 - BACK SECTION ASSEMBLY - MODEL 375
PAGE 11	- HEAD SECTION ASSEMBLY - MODEL 275
PAGE 12	- HEAD SECTION ASSEMBLY - MODEL 375
PAGE 13	- WRIST WREST ASSEMBLIES
PAGE 14	- PREVENTIVE MAINTENANCE GUIDE - OPTIONAL ACCESSORIES



IMPORTANT

REPORT ANY SHIPPING DAMAGE IMMEDIATELY !!
INFORM SHIPPER OF ANY DAMAGES — LEAVE
CARTON INTACT. LEAVE THE CHAIR IN THE
RECEIVING AREA UNTIL INSPECTION IS COMPLETE.

UNPACKING INSTRUCTIONS

Your A.P.C. series chair has been carefully packaged at our manufacturing plant to insure safe shipment to your medical facility. There are several procedures you must follow to put your new chair into service. These procedures only take a few minutes to complete and are required to insure proper operation of the chair.

-
- STEP 1 - CUT THE TWO METAL BANDS AROUND THE SHIPPING CARTON.
 - STEP 2 - REMOVE THE TOP HALF OF THE CARTON AND CUT ONE SIDE OF THE BOTTOM HALF.
 - STEP 3 - REMOVE THE CHAIR FROM THE CARTON.
 - STEP 4 - CHECK TO SEE IF ALL FEATURES OF THE CHAIR WORK PROPERLY. IF ALL FEATURES WORK GO TO STEP 5. IF ANY OF THE FEATURES OF THE CHAIR DO NOT WORK PROPERLY CALL HAUSTED FOR SERVICE AT: 1-800-428-7833 OR IN OHIO (216) 723-3271.
 - STEP 5 - CLEAN YOUR NEW CHAIR USING A MILD DETERGENT TO REMOVE ANY DIRT ACCUMULATED DURING SHIPMENT AND PLACE THE CHAIR INTO SERVICE.

HAUSTED[®] ASC - GENERAL INSTRUCTIONS & CHAIR SPECIFICATIONS

CHAIR SPECIFICATIONS

THE A.S.C. HAS BEEN DESIGNED TO MEET BOTH THE CLINICAL NEEDS FOR MINOR HEAD AREA TREATMENT AND PATIENT'S NEED FOR COMFORT.

PATIENT SURFACE SIZE		
CHAIR POSITION	STRETCHER POSITION	
21"	21"	WIDTH
39.5"	75"	LENGTH

UNIT HEIGHT		
CHAIR POSITION	STRETCHER POSITION	
57"	22.5"	MIN.
64.5"	30"	MAX.

A FEW FEATURES OF THE A.S.C. ARE:

- *SURGEON OPERATED HEADREST
- *HYDRAULIC HEIGHT ADJUSTMENT
- *TEN DEGREE TRENDLENBURG
- *4 WHEEL LOCK CASTER BRAKING
- *STEERING LOCK CASTER
- *PNEUMATIC BACKREST
- *INDEPENDENT LEG SECTION
- *ADJUSTABLE FOOT REST
- *BUILT IN STORAGE TRAY

PLUS:

- THE A.S.C. COMES STANDARD WITH:
 - REMOVABLE PRE-OP / POST-OP HEAD EXTENSION W/ PAD.
 - BUILT IN TELESCOPING SURGICAL SUPPORT LEGS.
 - THICK COMFORT CUSHIONS, FOREST GREEN IS STANDARD COLOR.
- * VARIETY OF DIFFERENT COLORED COMFORT CUSHIONS ARE AVAILABLE BY SPECIAL ORDER.

GENERAL INSTRUCTIONS

-AS A VALUED CUSTOMER AND USER OF HAUSTED'S UNIQUE A.S.C. CHAIR, THESE GENERAL INSTRUCTIONS ARE TO INSURE PROPER USE IS UNDERSTOOD AND FOLLOWED TO OBTAIN MAXIMUM SAFETY AND FUNCTION OF YOUR HAUSTED A.S.C. CHAIR.

-THE A.S.C. CHAIR'S UNEQUALED VERSATILITY FROM CHAIR POSITIONING TO FULL RECLINING IS DESIGNED TO PERMIT CONVENIENT PATIENT POSITIONING AND TREATMENT WHEN REQUIRED WITH A **MAXIMUM PATIENT WEIGHT OF 325 lbs.** THE A.S.C. **IS NOT** INTENDED TO REPLACE A STRETCHER OR GURNEY OF WHICH HAUSTED OFFERS A COMPLETE PRODUCT LINE FOR YOUR SPECIALTY NEEDS.

-THE A.S.C. HAS A "WARNING" LABEL AT THE **HEAD** AND **FOOT** END OF THE CHAIR FRAME, STATING: **DO NOT SIT ON END - AS TIPPING MAY OCCUR.**


-PATIENT **ENTRY, EGRESS** AND **TRANSFER** FROM THE A.S.C. CHAIR SHOULD ALWAYS BE FROM THE **CENTER SIDE RAIL LOCATIONS**, BY LOWERING OF THE SIDE RAILS AS SHOWN ABOVE. **AT NO TIME** SHOULD THE PATIENT BE PERMITTED TO ENTER OR EXIT FROM THE ENDS OF THE A.S.C. CHAIR WHEN THE CHAIR IS IN PARTIAL OR FULL RECLINE POSITION.

-**NOTE:** WHEN NOT IN USE, **DO NOT** LEAVE THE A.S.C. CHAIR IN A FULL OR PARTIAL RECLINE POSITION.

-**NOTE:** THE PATIENT SHOULD BE SECURED WITH RESTRAINT STRAPS WHEN LEFT UNATTENDED IN A FULL RECLINE POSITION.

-THESE GENERAL OPERATING INSTRUCTIONS ARE TO INSURE PATIENT SAFETY AND PERMITS YOU TO USE YOUR A.S.C. CHAIR TO TREAT, TRANSPORT AND TRANSFER YOUR PATIENT WITH UNEQUALED VERSATILITY AND SAFETY.

IF YOU HAVE ANY QUESTIONS, PLEASE CONTACT THE "HAUSTED CUSTOMER SERVICE DEPARTMENT" AT:
1-800-428-7833 OR (IN OHIO) 216-723-3271

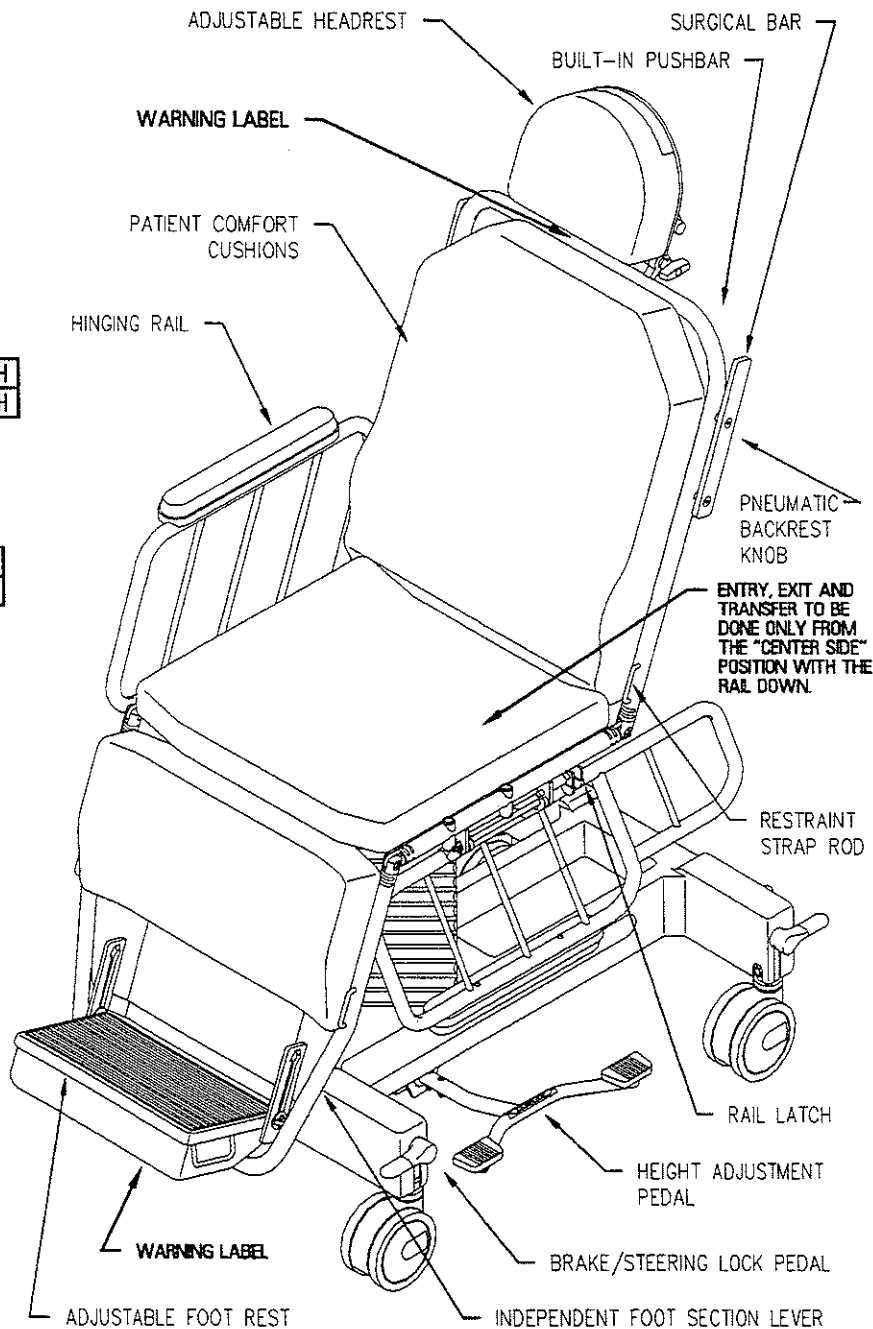

- OPERATING INSTRUCTIONS

BRAKING AND STEERING

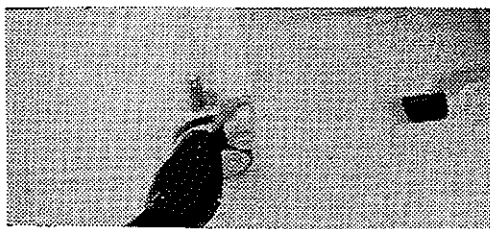


FIG. 1

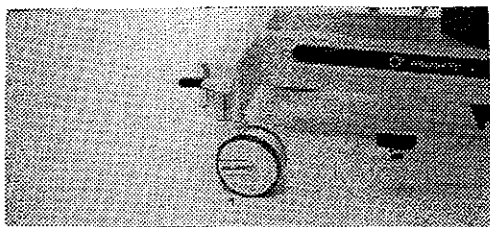


FIG. 2

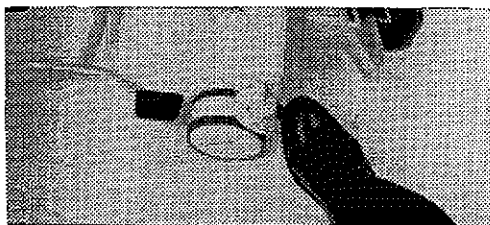


FIG. 3

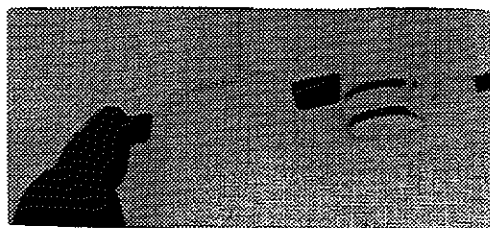


FIG. 4

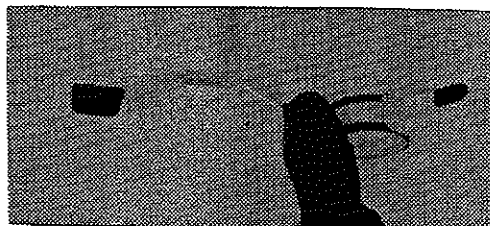


FIG. 5

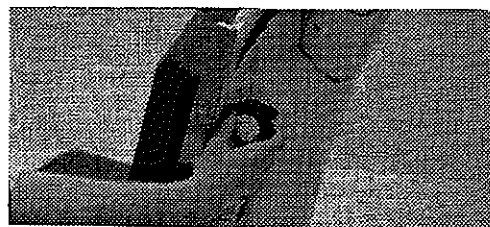


FIG. 6

APPLYING THE BRAKES

THE FOUR WHEEL CENTRAL BRAKING SYSTEM IS ACTIVATED BY DEPRESSING THE RED PEDAL AT EITHER THE LEFT OR RIGHT SIDE AT THE FOOT END OF THE UNIT (FIG. 1). TO FULLY ENGAGE THE BRAKES, THE PEDAL SHOULD BE PRESSED TO APPROXIMATELY 45°. ALL FOUR CASTER SHOULD BE LOCKED FROM SWIVELING AND ROTATING.

RELEASING THE BRAKES

DEPRESS THE BLACK PEDAL AT EITHER OF THE HEAD END CORNERS OF THE UNIT, UNTIL THE PEDAL IS IN A HORIZONTAL POSITION (FIG. 2). ALL FOUR WHEELS SHOULD ROTATE AND SWIVEL FREELY.

BRAKES SHOULD ALWAYS BE ACTIVATED WHEN UNIT IS NOT IN TRANSPORT

APPLYING THE STEERING LOCK

FROM THE HEAD END OF THE ASC CHAIR, DEPRESS THE BLACK PEDAL DOWNWARD INTO LOCKED POSITION (FIG. 3). PUSH THE ASC CHAIR FORWARD. ONE CASTER AT THE FOOT END WILL LOCK INTO A NON-SWIVEL MODE. THIS ALLOWS STRAIGHT STEERING BY THE ATTENDANT, WHEN PUSHING FROM THE HEAD END.

RELEASING THE STEERING LOCK

DEPRESS THE RED PEDAL AT EITHER OF THE FOOT END CORNERS OF THE UNIT, UNTIL THE PEDAL IS IN A HORIZONTAL POSITION (FIG. 2). ALL FOUR WHEELS SHOULD ROTATE AND SWIVEL FREELY.

TOP ADJUSTMENT

RAISING THE TOP

PRESS THE PUMP PEDAL (PEDAL NEAREST THE HEAD END) TO THE FLOOR (FIG. 4), THEN RELEASE THE PRESSURE ON THE PEDAL. REPEAT THIS PROCESS UNTIL DESIRED HEIGHT IS OBTAINED. USE SMOOTH STROKES ON PUMP PEDAL TO ENSURE PATIENT COMFORT.

LOWERING THE TOP

PRESS DOWN ON THE RELEASE PEDAL (PEDAL NEAREST THE FOOT END) (FIG. 5) UNTIL DESIRED HEIGHT IS OBTAINED.

RECLINING THE CHAIR INTO THE STRETCHER POSITION

ROTATE THE RECLINE KNOB, LOCATED ON EITHER SIDE OF THE CHAIR, TOWARD THE BACK SECTION (FIG. 6). PUSH BACK ON THE BACK UNTIL IT IS AT DESIRED POSITION, (FIG. 7) AND RELEASE THE KNOB. TO POSITION THE BACK INTO THE UPRIGHT POSITION, ROTATE THE RECLINE KNOB AND LIFT UP ON THE BACK SECTION.

THE A.S.C. SHOULD BE KEPT IN FULL UPRIGHT POSITION WHEN NOT IN USE

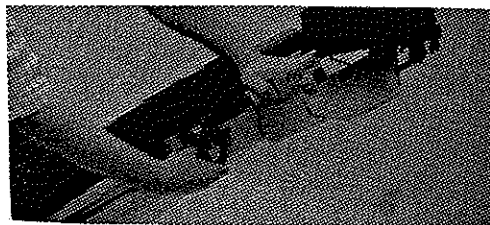


FIG. 7

- OPERATING INSTRUCTIONS

TOP ADJUSTMENT - CONTINUED

RECLINING THE CHAIR INTO TRENDLENBURG POSITION

ROTATE THE RECLINE KNOB, LOCATED ON EITHER SIDE OF THE CHAIR, AWAY FROM THE BACK SECTION (FIG. 8). PUSH BACK ON THE BACK SECTION UNTIL IT IS IN THE TRENDLENBURG POSITION (FIG. 9). AND RELEASE THE KNOB. TO REMOVE FROM THE TRENDLENBURG POSITION, ROTATE THE KNOB AWAY FROM THE CHAIR WHILE LIFTING UP ON THE BACK SECTION.

THE A.S.C. SHOULD BE KEPT IN FULL UPRIGHT POSITION WHEN NOT IN USE

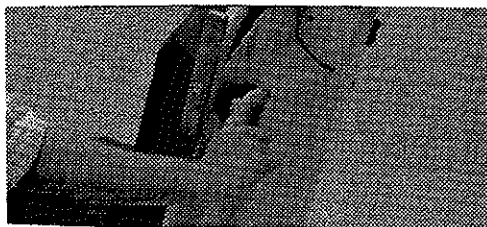


FIG. 8

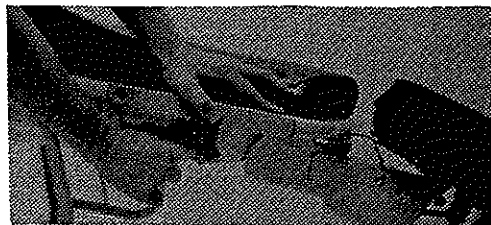


FIG. 9

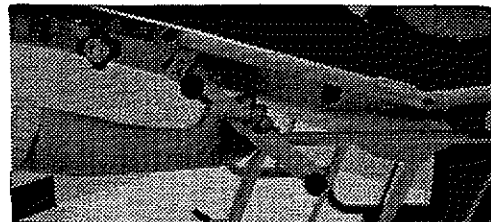


FIG. 10

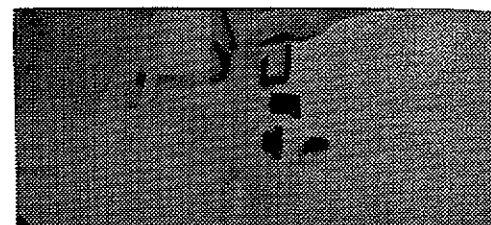


FIG. 11

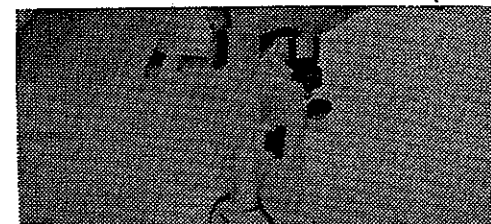


FIG. 12



FIG. 13

USING THE SURGICAL LEGS

ADJUSTING THE TOP

ADJUST THE TOP INCLINE TO DESIRED POSITION AND ADJUST HEIGHT TO DESIRED POSITION. RELEASE ONE OF THE SURGICAL LEGS FROM THE BACK OF THE CHAIR (FIG. 10) BY PULLING IT AWAY FROM THE BACK SECTION. EXTEND THE LEG TOWARD THE LEG SOCKET LOCATED ON THE BACK OF THE BASE ASSEMBLY (FIG. 11). IF THE LEG EXTENDS OVER THE SOCKET (**DO NOT FORCE THE LEG TO EXTEND BEYOND ITS BUILT IN SPRING**) (FIG. 12) THEN THE TOP IS ADJUSTED CORRECTLY, GO TO "ADJUSTING SURGICAL LEGS". IF THE LEG DOES NOT REACH THE SOCKET, (FIG. 11) YOU NEED TO READJUST THE TOP. YOU NEED TO LOWER THE HEIGHT OR DECREASE THE INCLINE UNTIL THE LEG REACHES THE SOCKET.

ADJUSTING SURGICAL LEGS

RELEASE THE SURGICAL LEGS FROM THE BACK OF THE CHAIR (FIG. 10) BY PULLING IT AWAY FROM THE BACK SECTION. EXTEND THE LEGS TOWARD THE LEG SOCKET (FIG. 11). (IF THE LEG DOES NOT REACH, THEN GO TO "ADJUSTING THE TOP" IN THIS SECTION) PLACE THE BOTTOM OF THE LEG OVER THE SOCKET (FIG. 12). TIGHTEN THE LOWER KNOB (FIG. 13) TO FIX THE LEG ONTO THE SOCKET. THEN TIGHTEN THE UPPER KNOB (FIG. 14) TO FIX KEEP THE LEG FROM TELESCOPING. (**ONCE THE LEGS ARE LOCKED INTO PLACE DO NOT TRY TO ADJUST THE TOP**)

DISENGAGING SURGICAL LEGS

LOOSEN THE LOWER KNOB FIRST (FIG. 13). SUPPORT THE LOWER LEG SECTION (FIG. 15) AND LOOSEN THE UPPER KNOB (FIG. 14). ALLOW THE LEG TO RETRACT AND PRESS BACK INTO THE STORAGE CLIP.

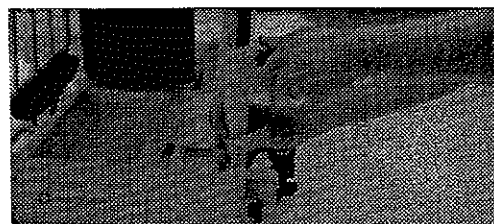


FIG. 14



FIG. 15

- OPERATING INSTRUCTIONS

INDEPENDENT LEG SECTION

INDEPENDENT OPERATION

HOLD AND SUPPORT LEG SECTION BY THE BOTTOM OF THE FRAME TUBE (FIG. 16). PULL OUT ON THE END OF THE RED LEVER (FIG. 16). ROTATE THE LEVER CLOCKWISE WHILE POSITIONING THE LEG SECTION (FIG. 17). ONCE DESIRED POSITION OF THE LEG SECTION IS REACHED, RELEASE THE LEVER AND THE LEG SECTION WILL BE RIDGED (FIG. 18). WHEN THE LEG SECTION IS INDEPENDENT, IT CAN BE REPOSITIONED AT ANY TIME BY TURNING THE LEVER CLOCKWISE (FIG. 17).

ALWAYS SUPPORT THE LEG SECTION WHEN TURNING THE RED LEVER

RETURNING TO CHAIR MODE

SUPPORT THE LEG SECTION, TURN THE LEVER COUNTER-CLOCKWISE WHILE PUSHING IN ON THE LEVER (FIG. 19). THE LEVER WILL SLIDE INWARD TO AN INACTIVE POSITION. ONCE THE LEVER IS INACTIVE, ADJUST THE LEG SECTION UNTIL IT IS IN LINE WITH THE BACK SECTION (FIG. 20). IT WILL CLICK AS IT RE-ENGAGES WITH THE BACK SECTION.

ADJUSTABLE FOOTREST

REPOSITIONING THE FOOTREST

THE FOOTREST HAS THREE POSITIONS; RETRACTED, LOWER AND THE UPPER POSITIONS. WITH THE PAN IN THE RETRACTED POSITION PULL OUT, ON THE HANDLES ON EITHER SIDE OF PAN (FIG. 21). THE FOOTREST WILL DROP INTO THE LOWER POSITION (FIG. 22). TO PLACE THE FOOTREST INTO THE UPPER POSITION, GRASP THE HANDLES ON BOTH SIDES OF THE PAN. TILT THE PAN UP WHILE LIFTING (FIG. 23). ONCE THE PAN IS UP ALL THE WAY, TILT THE PAN OUT UNTIL IT LOCKS INTO POSITION (FIG. 24). TO RETURN THE FOOTREST INTO THE RETRACTED POSITION, TILT THE PAN UP WHILE LETTING IT SLIDE DOWN, IF NEEDED, (FIG. 23), INTO THE LOWER POSITION (FIG. 21), CONTINUE TILTING THE FOOTREST UNTIL IT IS IN THE RETRACTED POSITION.

WARNING: NEVER STAND ON THE FOOTREST

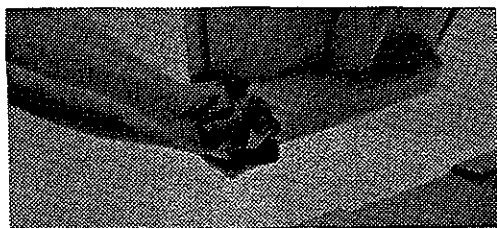


FIG. 16

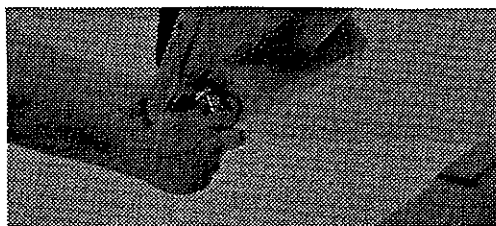


FIG. 17

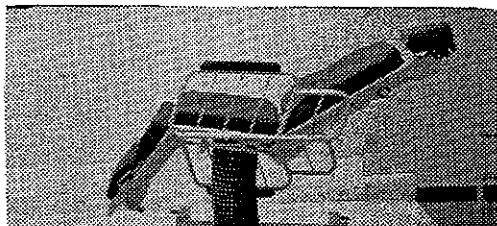


FIG. 18

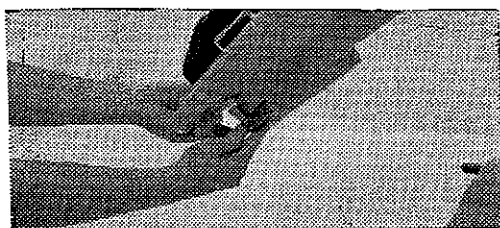


FIG. 19

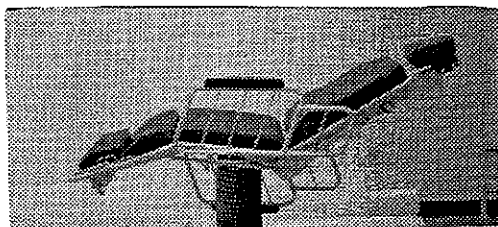


FIG. 20

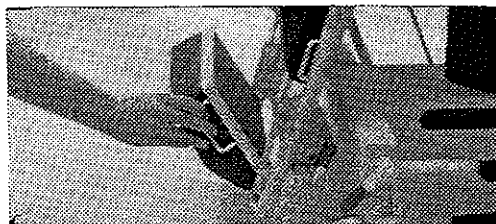


FIG. 21

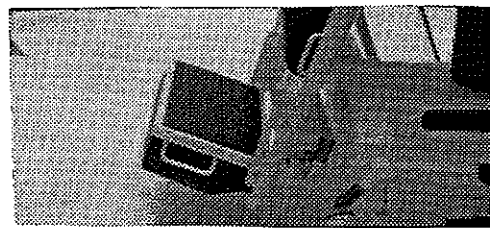


FIG. 22

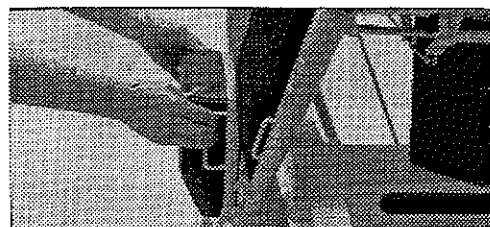


FIG. 23

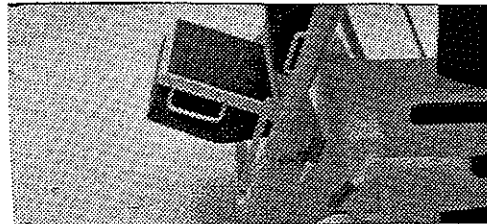


FIG. 24

- OPERATING INSTRUCTIONS

HINGING RAILS

REPOSITIONING THE RAIL

THE APC RAIL HAS TWO POSITIONS, UP AND DOWN, BOTH POSITIONS LOCK THE RAIL INTO PLACE. TO PUT THE RAIL DOWN, WHILE SUPPORTING THE RAIL (FIG. 25), SQUEEZE THE TWO RED TABS UNDER THE SEAT (FIG. 26). ALLOW THE RAIL TO HINGE ALL THE WAY DOWN UNTIL IT LOCKS INTO PLACE (FIG. 27). TO RAISE THE RAIL, SQUEEZE THE TWO RED TABS UNDER THE SEAT (FIG. 26). LIFT THE RAIL ALL THE WAY UP UNTIL IT LOCKS INTO PLACE.

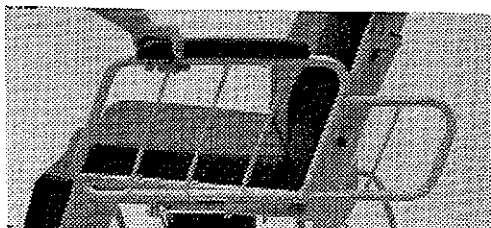


FIG. 25

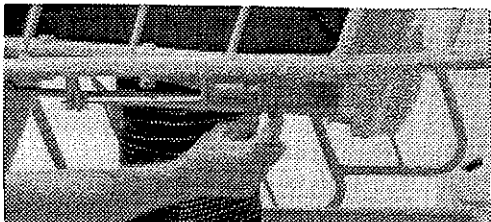


FIG. 26

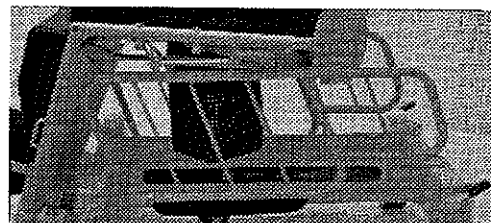


FIG. 27



FIG. 28



FIG. 29

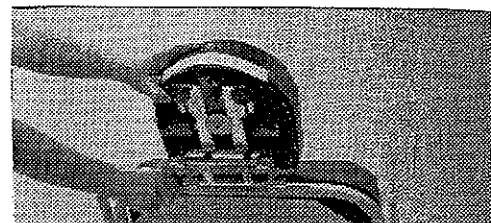


FIG. 30

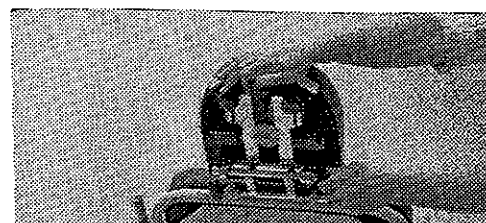


FIG. 32



FIG. 33

A.S.C. HEADREST

PRE-OP / POST-OP HEAD EXTENSION

TO REMOVE THE HEAD EXTENSION, LOOSEN THE "T" KNOB UNDER THE HEADREST (FIG. 28). WITH THE KNOB LOOSE PULL OUT ON THE HEAD EXTENSION. TO REPLACE THE EXTENSION, SLIDE THE EXTENSION OVER THE HEAD PIECE USING GUIDES. PUSH THE EXTENSION AGAINST THE HEADREST AND TIGHTEN "T" KNOB (FIG. 28).

ADJUSTING THE ASC 275 HEADREST

REMOVE HEAD EXTENSION SEE "PRE-OP/POST-OP HEAD EXTENSION" INSTRUCTIONS. GRASP LARGE T-KNOB ON BOTTOM OF HEADREST AND GRASP SIDE OF HEAD REST (FIG. 29). PULL T-KNOB TOWARD THE BACK OF HEADREST. WITH THE T-KNOB PULLED BACK HEADREST PIVOT HEADREST TO DESIRED POSITION THEN RELEASE THE KNOB.


ADJUSTING THE ASC 375 HEADREST

REMOVE HEAD EXTENSION, SEE "PRE-OP/POST-OP HEAD EXTENSION" INSTRUCTIONS. LOOSEN THE KNOB ON THE LEFT SIDE OF THE HEAD REST AND POSITION THE HEIGHT OF THE HEADREST (FIG. 30). ONCE THE DESIRED POSITION IS REACHED, TIGHTEN THE KNOB AND TEST THE HEADREST TO ENSURE RATCHET TEETH ARE LOCKED. LOOSEN THE KNOB ON THE RIGHT SIDE OF THE HEADREST AND POSITION THE ANGLE OF THE HEADREST (FIG. 31). ONCE THE DESIRED ANGLE IS REACHED, TIGHTEN THE KNOB AND TEST THE HEADREST TO ENSURE RATCHET TEETH ARE LOCKED. LOOSEN THE TWO KNOBS ON THE BOTTOM OF THE HEADREST AND ADJUST HEADREST LENGTH (FIG 32), TIGHTEN KNOBS.

BEFORE PLACING PATIENTS HEAD ON HEADREST, APPLY PRESSURE TO THE *HEADREST TO ENSURE LOCKING OF THE RATCHET TEETH

HAUSTED® ASC

MODEL 275 AND 375 BASE ASSEMBLY

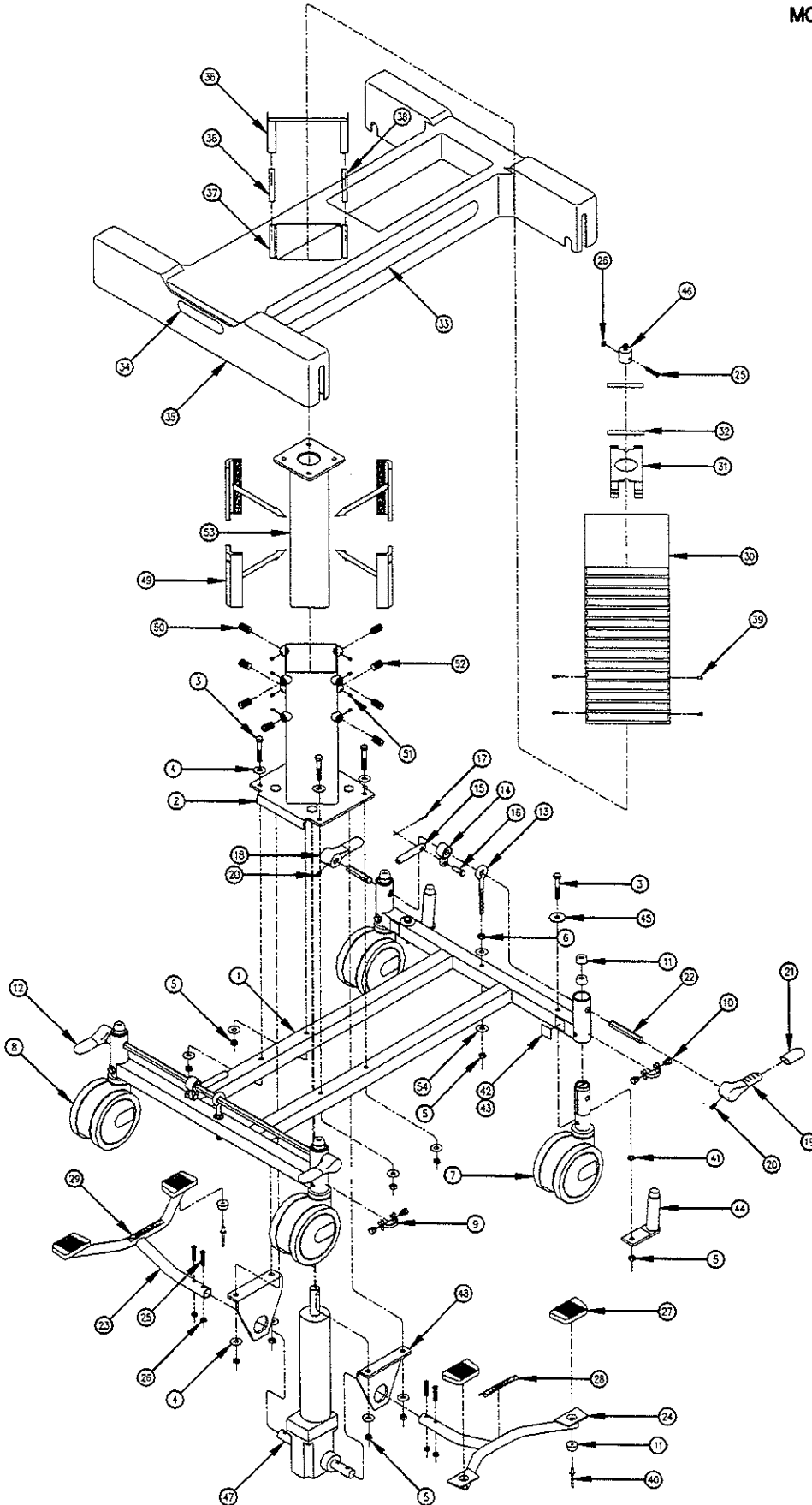


54	031080	5/16 FLAT WASHER	4
53	068867	INSIDE COLUMN ASSEMBLY	1
52	068199	COLUMN GLIDE SCREW W/O POINT	4
51	052041	# 10-24 x 1/4 CUP POINT SET SCR	8
50	061116	COLUMN GLIDE SCREW	4
49	068178	SHORT GLIDE SUPPORT	4
48	061122	HYDRAULIC JACK SUPPORT HANGER	2
47	061122	MEDBACK	1
46	067354	APC SEAT MOUNTING BOSS	1
45	031359	3/8 I.D. x 1-1/4 O.D. FLAT WASHER	2
44	068306	SURGICAL LEG STOP WELDMENT	2
43	061155	VELCRO HOOK FOR COVER	4
42	061156	VELCRO PILE FOR COVER	4
41	031092	5/16 INTERNAL TOOTH LOCK WASHER	2
40	000794	POP RIVET	2
39	065129	BELLOWS CUP	4
38	061166	VELCRO PILE 5"	2
37	066167	VELCRO HOOK 4"	2
36	061287	BELLOWS MOUNTING ANGLE ASSY	1
35	065399	BASE COVER	1
34	061192	HAUSTED PRODUCT LABEL	2
33	061193	APC CHAIR LABEL	2
32	061168	TRIM-LOC 3-1/2"	1
31	061379	BELLOWS TOP PLATE	1
30	061261	APC BELLOWS	1
29	061229	PUMP PEDAL LABEL - LEFT	1
28	061228	PUMP PEDAL LABEL - RIGHT	1
27	061164	PUMP PEDAL PAD	4
26	031055	1/4-20 LOCK NUT	5
25	062128	1/4-20 x 1-1/4 BUTTON HD SOCKET SCR	5
24	068045	PUMP PEDAL - RIGHT	1
23	068046	PUMP PEDAL - LEFT	1
22	061127	HEX BRAKE ROD	2
21	061184	BRAKE PEDAL COVER - BLACK	2
20	031101	1/4-20 x 1/2 SET SCR	2
19	061183	BRAKE PEDAL - LEFT	2
18	061182	BRAKE PEDAL - RIGHT	2
17	031466	3/32 x 1/2 COTTER PIN	2
16	052419	LINK PIN	2
15	065366	CONNECTING LINK	1
14	052684	CAM LINK - SHORT	2
13	062688	EYE BOLT	2
12	061183	BRAKE PEDAL CAP - RED	2
11	031156	RUBBER GROMMET	10
10	062685	M8 x 10 HEX HD CAP SCR	8
9	062659	CASTER BOLT LOCK PLATE	4
8	065294	TWIN WHEEL STEERING LOCK	1
7	065293	TWIN WHEEL TOTAL LOCK CASTER	1
6	031058	5/16-18 FREE THREAD HEX NUT	2
5	031057	5/16-18 LOCK NUT	12
4	031076	1/4 SAE FLAT WASHER	12
3	031334	5/16-18 x 2 HEX HD CAP SCR	6
2	061117	OUTSIDE COLUMN ASSEMBLY	1
1	065398	APC BASE FRAME ASSY	1

ITEM	PART NO	DESCRIPTION	QTY
------	---------	-------------	-----

HAUSTED® ASC


MODEL 275 AND 375 SEAT ASSEMBLY

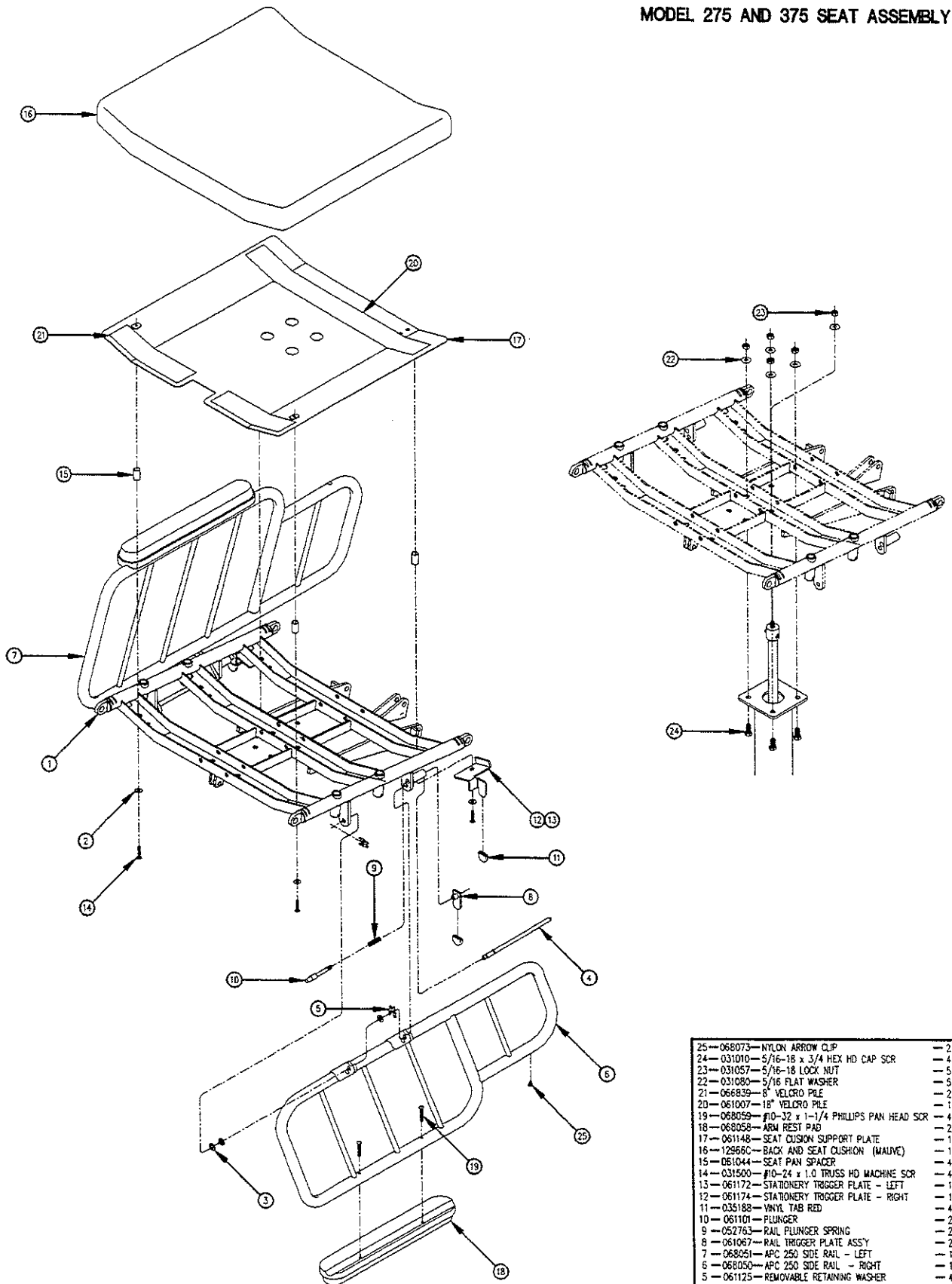

25	068073	NYLON ARROW CLIP	2
24	031010	5/16-18 x 3/4 HEX HD CAP SCR	4
23	031057	5/16-18 LOCK NUT	5
22	031080	5/16 FLAT WASHER	5
21	066839	8" VELCRO PILE	2
20	061007	18" VELCRO PILE	1
19	068059	#10-32 x 1-1/4 PHILLIPS PAN HEAD SCR	4
18	068058	ARM REST PAD	1
17	061148	SEAT CUSHION SUPPORT PLATE	1
16	129666	BACK AND SEAT CUSHION (MAINE)	1
15	061044	SEAT PAN SPACER	4
14	031500	#10-24 x 1.0 TRUSS HD MACHINE SCR	4
13	061172	STATIONERY TRIGGER PLATE - LEFT	1
12	061174	STATIONERY TRIGGER PLATE - RIGHT	1
11	035188	VINYL TAB RED	4
10	061101	PLUNGER	2
9	052763	RAIL PLUNGER SPRING	2
8	061067	RAIL TRIGGER PLATE ASS'Y	2
7	068051	APC 250 SIDE RAIL - LEFT	1
6	068050	APC 250 SIDE RAIL - RIGHT	1
5	061125	REMOVABLE RETAINING WASHER	4
4	061092	RAIL PIVOT ROD	2
3	035115	3/4 OD x .395 ID NYLON WASHER	6
2	031364	1/4 SAE FLAT WASHER	4
1	068039	APC 5" SEAT SECTION WELDMENT	1

HAUSTED® ASC


MODEL 275 AND 375 FOOT SECTION

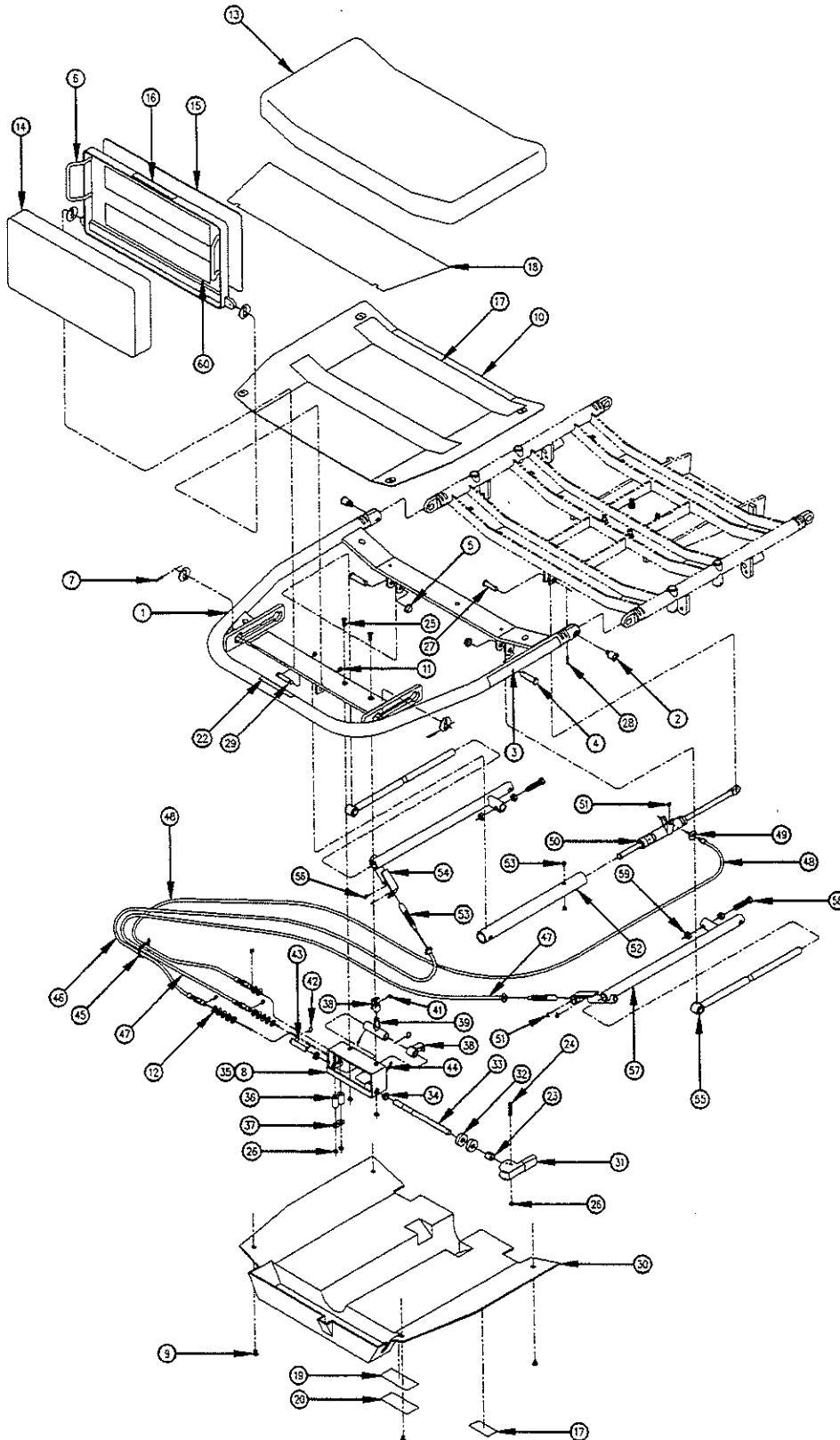

60	061051	2.0 x 11.0 VELCRO PILE	2
59	031057	5/16-18 LOCK NUT	4
58	061350	5/16-18 x 1-1/2 HEX HD CAP SCR	2
57	061325	PLUNGER TUBE ASSEMBLY - RIGHT	1
56	061324	3/8 INTERNAL RETAINING RING	2
55	061328	PLUNGER ROD ASSEMBLY	2
54	061320	PLUNGER TUBE ASSEMBLY - LEFT	1
53	031102	1/4-20 x 1/4 SET SCREW - KNURLED	2
52	061331	LEG EXTENSION TUBE	1
51	062641	RETAINING RING 5103-25	3
50	060105	MECHLOK - MECHANICAL LOCK	1
49	031364	1/4 SAE FLAT WASHER	1
48	061310	A.P.C. LEG CONTROL CABLE	1
47	061311	RIGHT PLUNGER CABLE	1
46	061312	LEFT PLUNGER CABLE	1
45	000849	CABLE TIE	1
44	031111	1/8 x 5/8 SPRING PIN	2
43	061339	SUPPORT SHAFT	1
42	031137	RETAINING RING 5103-37	2
41	051344	RETAINING RING 5103-31	4
40	061342	CABLE PULL ASSEMBLY	1
39	061346	CONTROL PIVOT SUB-ASSEMBLY	1
38	061348	PULL TUBE ASSEMBLY	1
37	031074	3/16 SAE WASHER	2
36	061337	ROLLER	2
35	061333	CONTROL BOX BRACKET ASSEMBLY	1
34	061336	FLANGE NYLINER 6L2-FF	2
33	061338	CONTROL SHAFT	1
32	065411	ACTIVATING LEVER BEARING	2
31	065409	ACTIVATING LEVER	1
30	061306	LEG SECTION COVER	1
29	061351	LEG EXTENSION TUBE PIN	1
28	031466	9/32 x 1/2 COTTER PIN	1
27	060018	5/16 x 1.0 LINK PIN	1
26	031421	#10-24 LOCK NUT	5
25	031432	#10-24 x 5/8 FLAT SOCKET HD SCR	2
24	031500	#10-24 x 1.0 TRUSS HD MACH SCR	1
23	065412	ACTIVATING LEVER SPACER	1
22	061188	HEAD AND FOOT WARNING LABEL	2
20	031458	MYLAR COVER FOR IDENTIFICATION TAG	1
19	031457	SERIAL IDENTIFICATION TAG	1
18	068068	FOOT SECTION VINYL COVER	1
17	061007	2.0 x 18.0 VELCRO PILE	2
16	061141	FOOT REST WARNING LABEL	1
15	068113	MAT FOR FOOT REST	1
14	12968E	FOOT REST CUSHION (MAUVE)	1
13	12966F	LEG SECTION CUSHION (MAUVE)	1
12	031081	SHIM WASHER	12
11	031137	RETAINING RING	2
10	060072	LEG PAD SUPPORT PLATE	1
9	061196	#10-24 x 3/8 TRUSS HD MACHINE SCR	4
8	061302	FOOT END CONTROL ASSY	1
7	061458	1/8 x 3/4 SPRING PIN	2
6	068060	FOOT REST ASSY	1
5	030287	3/8 PAL NUT	2
4	011318	3/8 x 1-7/16 CLEVIS PIN	2
3	061358	OPERATING LABEL	1
2	031036	5/16-24 x 15/16 PHILLIPS OVAL HD SCR	2
1	068035	FOOT SECTION ASSY	1

HAUSTED® ASC

MODEL 275 AND 375 BACK SECTION

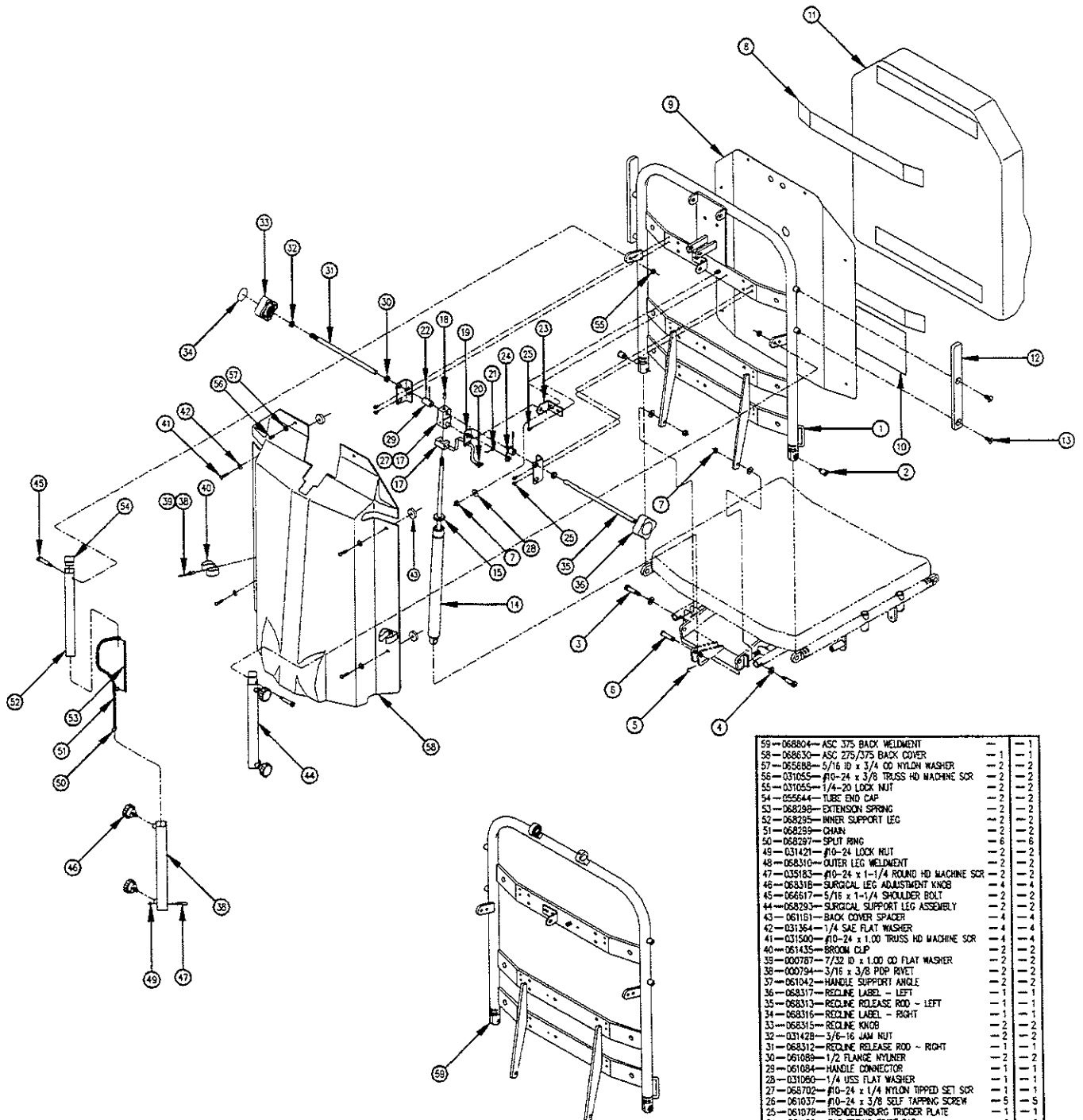


59	068804	ASC 375 BACK WELDMENT	1	1
58	068630	ASC 275/375 BACK COVER	1	1
57	065888	5/16 ID x 3/4 OD NYLON WASHER	2	2
56	031055	#10-24 x 3/8 TRUSS HD MACHINE SCR	2	2
55	031055	1/4-20 LOCK NUT	2	2
54	055644	TUBE END CAP	2	2
53	068286	EXTENSION SPRING	2	2
52	068295	INNER SUPPORT LEG	2	2
51	068299	CHAIN	2	2
50	068297	SPLIT RING	2	2
49	031421	#10-24 LOCK NUT	2	2
48	068310	OUTER LEG WELDMENT	2	2
47	035183	#10-24 x 1-1/4 ROUND HD MACHINE SCR	2	2
46	068316	SURGICAL LEG ADJUSTMENT KNOB	4	4
45	066617	5/16 x 1-1/4 SHOULDER BOLT	2	2
44	068293	SURGICAL SUPPORT LEG ASSEMBLY	2	2
43	061181	BACK COVER SPACER	4	4
42	031364	1/4 SAE FLAT WASHER	4	4
41	031500	#10-24 x 1.00 TRUSS HD MACHINE SCR	4	4
40	061435	BROOM CLIP	2	2
39	000767	7/32 ID x 1.00 OD FLAT WASHER	2	2
38	000794	3/16 x 3/8 POP RIVET	2	2
37	061042	HANDLE SUPPORT ANGLE	2	2
36	068317	RECLINE LABEL - LEFT	1	1
35	068313	RECLINE RELEASE ROD - LEFT	1	1
34	068316	RECLINE LABEL - RIGHT	1	1
33	068315	RECLINE KNOB	2	2
32	031428	3/6-16 JAM NUT	2	2
31	068312	RECLINE RELEASE ROD - RIGHT	1	1
30	061088	1/2 FLANGE NYLONER	2	2
29	061084	HANDLE CONNECTOR	1	1
28	031060	1/4 USS FLAT WASHER	1	1
27	068702	#10-24 x 1/4 NYLON TIPPED SET SCR	1	1
26	061037	#10-24 x 3/8 SELF TAPPING SCREW	5	5
25	061078	TRENDLENBURG TRIGGER PLATE	1	1
24	061152	GAS SPRING PIVOT EAR	1	1
23	031111	1/8 x 5/8 SPRING PIN	3	3
22	061063	TRENDLENBURG RELEASE SPRING	1	1
21	031014	#10-24 x 3/8 SOCKET HD CAP SCR	2	2
19	061099	TRENDLENBURG STOP PLATE WELDMENT	1	1
18	061153	ACTUATOR PIN	1	1
17	068701	APC ACTUATOR CASTING	1	1
16	061083	TRENDLENBURG STOP BLOCK	1	1
15	068900	GAS SPRING SPACER	1	1
14	067287	GAS SPRING	1	1
13	031034	5/16-18 x 3/4 FLAT SOCKET HD CAP SCR	4	4
12	068300	SURGICAL BAR	2	2
11	068626	ACS 275/375 BACK PAD	1	1
10	061191	APC BACK/SEAT VINYL COVER	1	1
9	068629	ACS 275/375 BACK SUPPORT PAN	1	1
8	061007	2.0 x 18.0 VELCRO PILE	2	2
7	031057	5/16-18 LOCK NUT	2	2
6	061296	GAS SPRING HINSE PIN	1	1
5	031137	RETAINING RING	1	1
4	031115	3/4 OD x .395 ID NYLON WASHER	4	4
3	061157	3/8 x 1.000 SHOULDER BOLT	2	2
2	031036	5/16-24 x 15/16 PHILLIPS OVAL HD SCR	2	2
1	068035	ACS 275 BACK SECTION WELDMENT	1	1

ITEM	PART NO	DESCRIPTION	QTY 275	QTY 375
------	---------	-------------	---------	---------

HAUSTED® ASC

MODEL 275 HEAD SECTION

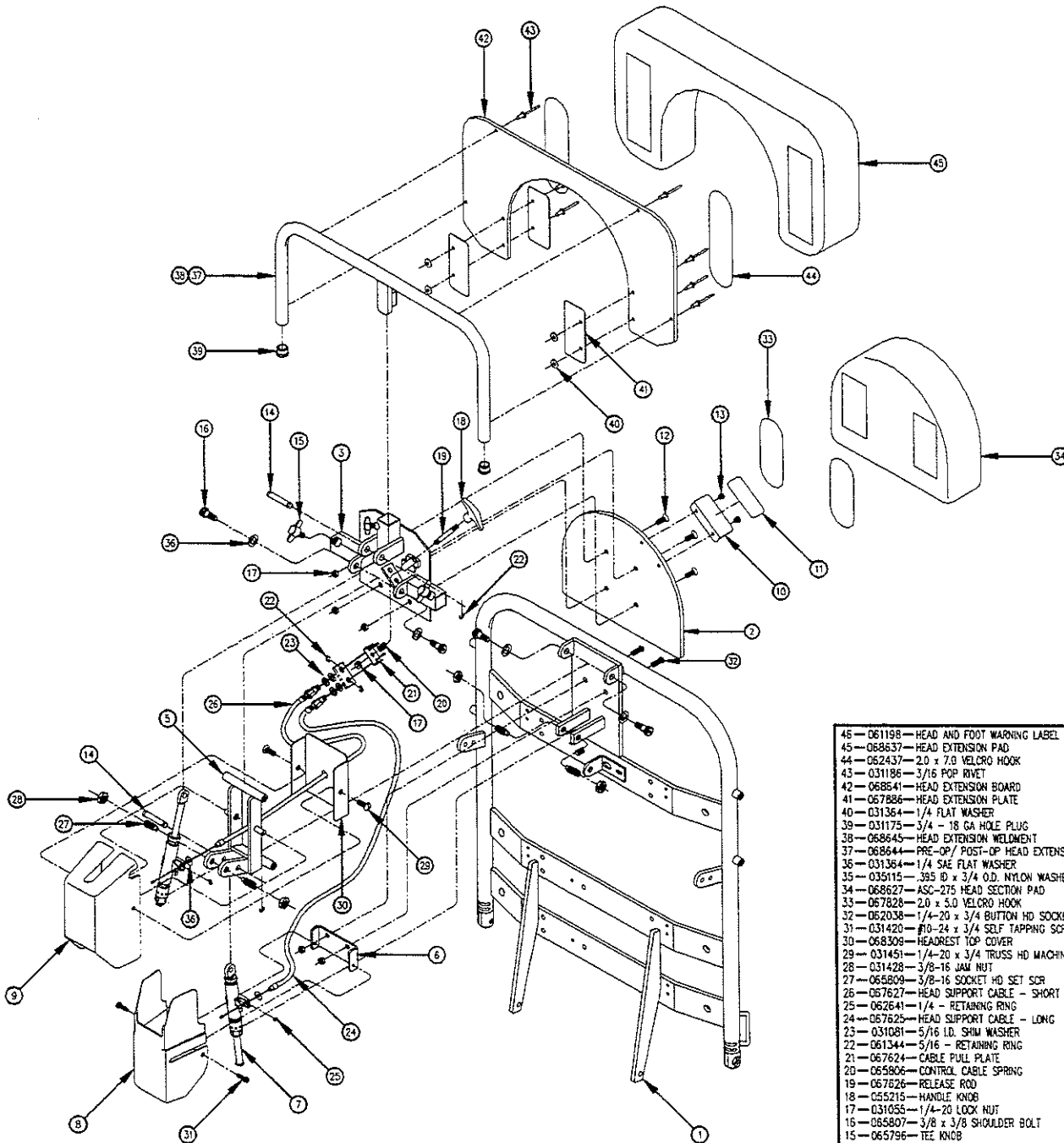


45	061198	HEAD AND FOOT WARNING LABEL (N.S.)	1
45	068637	HEAD EXTENSION PAD	1
44	062437	2.0 x 7.0 VELCRO HOOK	2
43	031186	3/16 POP RIVET	8
42	068541	HEAD EXTENSION BOARD	1
40	067886	HEAD EXTENSION PLATE	4
36	031364	1/4 FLAT WASHER	4
35	031175	5/4 - 18 GA HOLE PLUG	2
38	068645	HEAD EXTENSION WELDMENT	1
37	068644	PRE-OP/ POST-OP HEAD EXTENSION	1
36	031364	1/4 SAE FLAT WASHER	2
35	035115	.385 ID x 3/4 O.D. NYLON WASHER	4
34	068627	ASC-275 HEAD SECTION PAD	1
33	067828	2.0 x 5.0 VELCRO HOOK	2
32	062038	1/4-20 x 3/4 BUTTON HD SOCKET SCR	2
31	031420	#10-24 x 3/4 SELF TAPPING SCR	2
30	068308	HEADREST TOP COVER	1
29	031451	1/4-20 x 3/4 TRUSS HD MACHINE SCR	2
28	031428	3/8-16 JAM NUT	4
27	065809	3/8-16 SOCKET HD SET SCR	4
26	067627	HEAD SUPPORT CABLE - SHORT	1
25	062641	1/4 - RETAINING RING	2
24	067625	HEAD SUPPORT CABLE - LONG	1
23	031081	5/16 I.D. SHIM WASHER	4
22	061344	5/16 - RETAINING RING	4
21	067624	CABLE PULL PLATE	1
20	065806	CONTROL CABLE SPRING	1
19	067626	RELEASE ROD	1
18	055215	HANDLE KNOB	1
17	031053	1/4-20 LOCK NUT	6
16	065807	3/8 x 3/8 SHOULDER BOLT	4
15	065796	TEE KNOB	3
14	065808	HEAD SUPPORT PIN	2
13	031326	#10-24 x 1/4 BUTTON HD SOCKET SCR	2
12	031320	1/4-20 x 3/4 FLAT SOCKET HD CAP SCR	3
11	065802	OPERATING LABEL	1
10	065813	OPERATING LABEL ANGLE	1
9	057618	HEAD FRAME COVER	1
8	068314	TOP COVER	1
7	067622	MECHLOCK	2
6	067620	COVER MOUNTING PLATE	1
5	067614	HEAD SUPPORT FRAME	1
4	031324	1/4-20 x 1/2 FLAT SOCKET HD CAP SCR	1
3	068531	SUPPORT PLATE WELDMENT	1
2	058632	HEAD BOARD	1
1	068624	ASC275 BACK SECTION ASSY	1

ITEM	PART NO	DESCRIPTION	QTY
------	---------	-------------	-----

HAUSTED® ASC


MODEL 375 HEAD SECTION

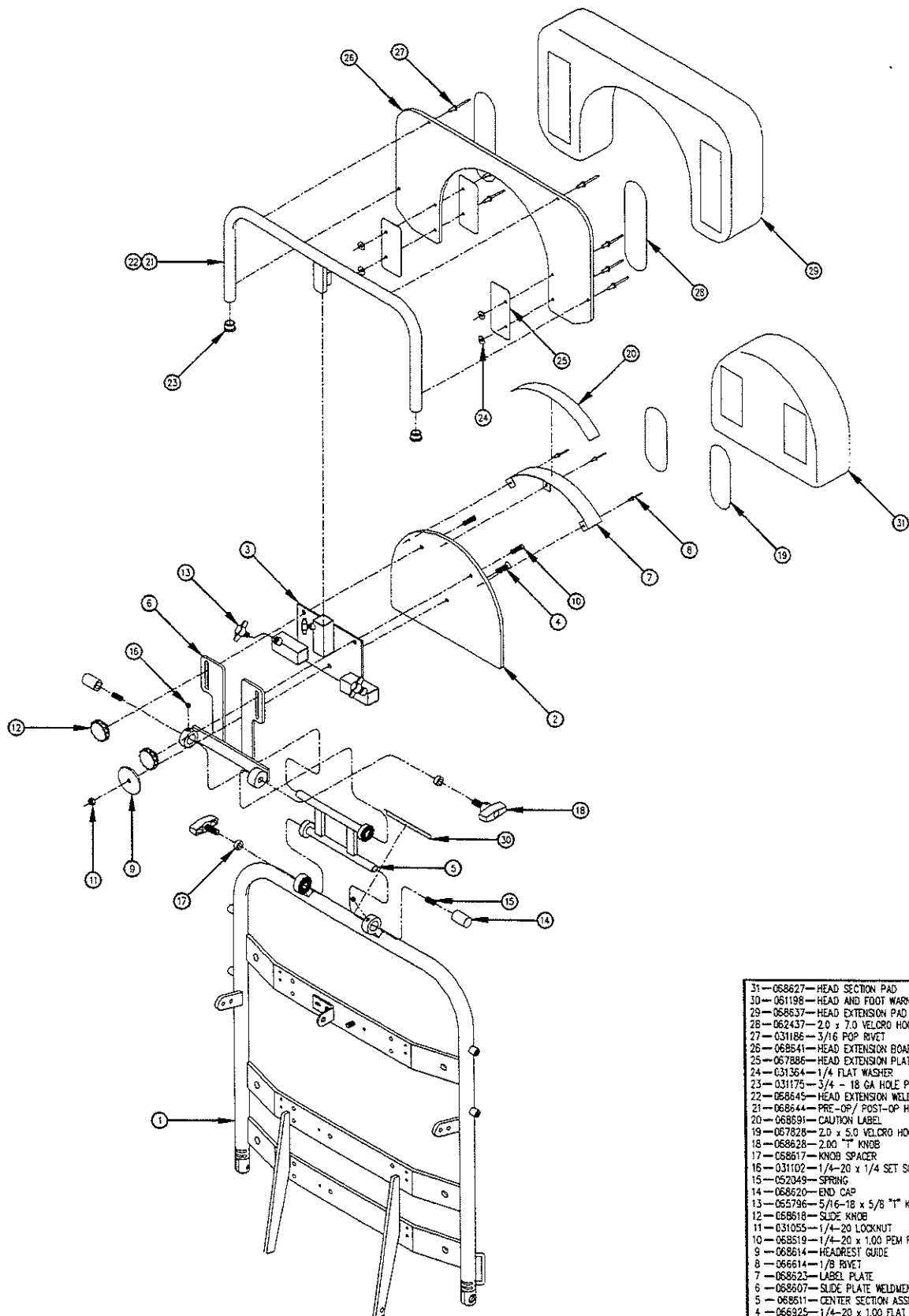

31	068627	HEAD SECTION PAD	1
30	061198	HEAD AND FOOT WARNING LABEL	1
29	068637	HEAD EXTENSION PAD	1
28	062437	2.0 x 7.0 VELCRO HOOK	2
27	031186	3/16 POP RIVET	3
26	068641	HEAD EXTENSION BOARD	1
25	067886	HEAD EXTENSION PLATE	4
24	031364	1/4 FLAT WASHER	4
23	031175	3/4 - 18 GA HOLE PLUG	2
22	068645	HEAD EXTENSION WELDMENT	1
21	068644	PRE-OP/ POST-OP HEAD EXTENSION	1
20	068691	CAUTION LABEL	1
19	067828	2.0 x 5.0 VELCRO HOOK	2
18	068628	2.00 "I" KNOB	2
17	068617	KNOB SPACER	2
16	031102	1/4-20 x 1/4 SET SCREW	2
15	052049	SPRING	2
14	068620	END CAP	2
13	065796	5/16-18 x 5/8 "I" KNOB ASSEMBLY	3
12	068616	SLIDE KNOB	2
11	031055	1/4-20 LOCKNUT	1
10	068519	1/4-20 x 1.00 PEM FASTENER	2
9	068614	HEADREST GUIDE	1
8	066614	1/8 RIVET	3
7	068623	LABEL PLATE	1
6	068607	SLIDE PLATE WELDMENT	1
5	068611	CENTER SECTION ASSEMBLY	1
4	066925	1/4-20 x 1.00 FLAT SOCKET HD CAP SCR	1
3	068615	SLIDE SUPPORT PLATE WELDMENT	1
2	068622	ASC375 HEAD BOARD	1
1	068604	ASC375 BACK SECTION ASSY	1
ITEM	PART NO	DESCRIPTION	QTY


HAUSTED® ASC

WRIST WREST ASSEMBLIES

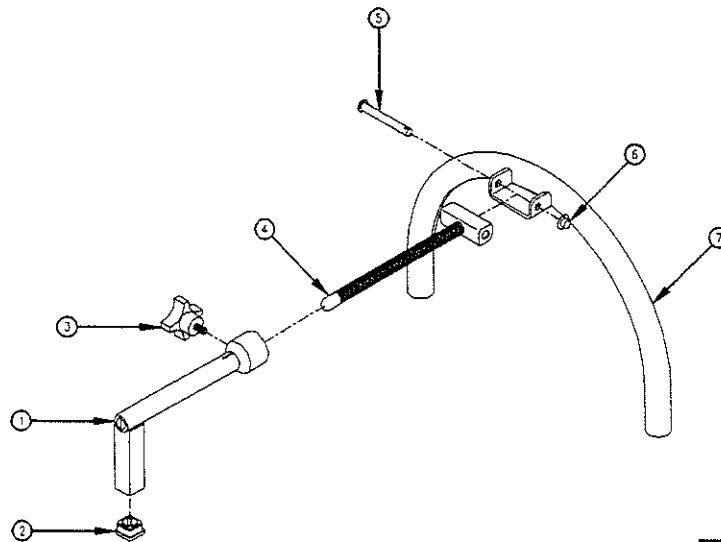

ACCESSORY NO: SA0746

7	065959	WRIST REST SUB-ASSEMBLY	1
6	031444	PAL CAP	1
5	056954	CLEVIS PIN	1
4	066740	HINGE BAR ASSEMBLY	1
3	066738	KNOB	1
2	067849	SQUARE INSERT	1
1	067736	WRIST REST SUPPORT ASSEMBLY	1
ITEM	PART NO	DESCRIPTION	QTY

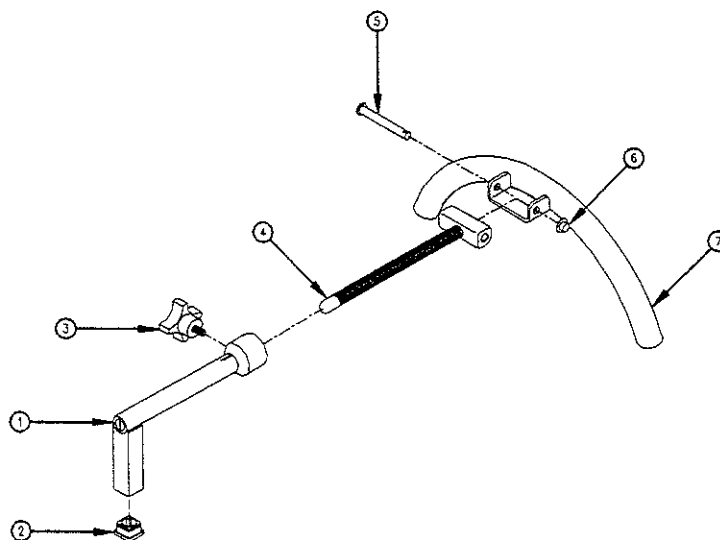

ACCESSORY NO: SA0729

7	067634	WRIST REST SUB-ASSEMBLY	1
6	031444	PAL CAP	1
5	056954	CLEVIS PIN	1
4	066740	HINGE BAR ASSEMBLY	1
3	066738	KNOB	1
2	067849	SQUARE INSERT	1
1	067736	WRIST REST SUPPORT ASSEMBLY	1
ITEM	PART NO	DESCRIPTION	QTY



RECOMMENDED PREVENATIVE MAINTENANCE

PROCEDURE	SCHEDULE	MATERIAL
LUBRICATE ALL MOVING AND SLIDING PARTS AND HINGE POINTS.	6 MONTHS	LUBRICATING OIL, LIGHT DUTY GREASE, WAX STICK LUBRICANT, OR NEVER-SEEZ LUBRICANT.
INSPECT ALL FASTENERS TO INSURE PROPER FIT, POSITION, AND TIGHTNESS (INCLUDING NUTS, BOLTS, ETC...)	3 MONTHS	PROPER SIZE WRENCH AND SCREWDRIVER.
CLEANING MATTRESS AND PADS	PERIODICALLY	SOAP, WATER, AND MEDIUM STRENGTH CLEANERS. (DO NOT PRESSURE CLEAN)
INSPECT ALL SURFACES AND REMOVE ANY USER DEVELOPED SHARP OR BURRED AREAS. APPLY TOUCH-UP WHERE REQUIRED.	PERIODICALLY	METAL FILE, PROPER COLOR PAINT (SPECIFY COLOR WHEN ORDERING)
CLEAN EXTERNAL SURFACES TO REMOVE DIRT, LINT, BLOOD, AND ANY OTHER DEPOSITS.	PERIODICALLY	MEDIUM STRENGTH ANTI-BACTERIAL CLEANERS, STAINLESS STEEL CLEANER & POLISH

NOTE:

IF UNITS ARE PRESSURE WASHED, LUBRICATION OF ALL MOVING AND SLIDING PARTS SHOULD BE DONE AFTERWARDS.

(STEAM CLEANING IS NOT RECOMMENDED AND CAN VOID WARRANTY)

OPTIONAL ACCESSORIES

NUMBER	DESCRIPTION	NUMBER	DESCRIPTION
SA0100	SURGICAL RAIL	APC018	TELESCOPING I.V. ROD, 27" TO 50" 2 SECTION, WITH HOLDER, STAINLESS
SA0201	ACCESSORY RAIL CLAMP	APC026	HEEL STIRRUPS WITH MOUNTING ADAPTORS
SA0315	KNEE CRUTCH LEG HOLDER WITH MOUNTING ADAPTOR	APC027	KNEE CRUTCHES WITH MOUNTING ADAPTORS
SA0414	ARMBOARD WITH 1" PAD AND CLAMP	APC14D	OXYGEN TANK HOLDER
SA0502	SHOULDER ARTHROSCOPY SYSTEM	APCSFT	CUSTOM COLOR COMFORT CUSHIONS SPECIAL ORDER, FOREST GREEN STD
SA0548	DOUBLE PKG. ARMTRAPS 6/CS	061280	HEAVY DUTY 5" CASTERS
SA0563	SHOULDER ARTHROSCOPY SYSTEM - WEIGHTLESS	088167	PATIENT TRAY
SA0729	LATERAL ACCESS SURGEON WRIST WREST	128450	VERTICAL OXYGEN TANK HOLDER
SA0746	FULL SURGEON (OVER THE BROW) WRIST WREST	130080	PAPER ROLL HOLDER
APC012	RESTRAINT STRAPS WITH VELCRO		

FOR FURTHER ASSISTANCE IN PREVENATIVE MAINTENANCE OR TO ORDER ACCESSORIES PLEASE CONTACT

HAUSTED

927 LAKE ROAD, P.O. BOX 710 MEDINA, OHIO 44258
(216) 723-3217 1-800-428-7833